


Object Configuration

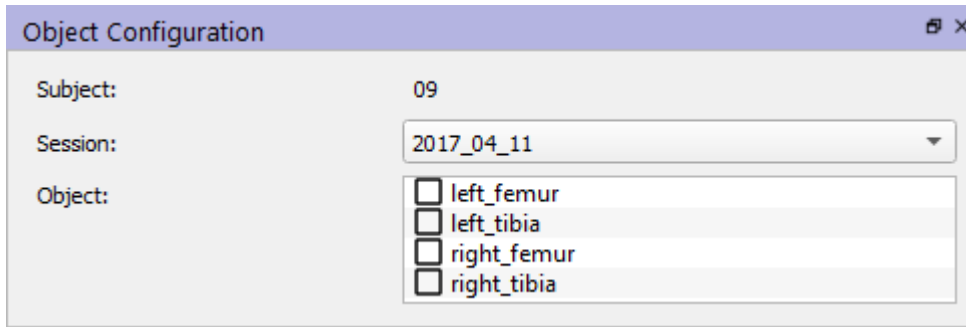

The screenshot shows a window titled "Object Configuration". It contains three main sections: "Subject:" with the value "09", "Session:" with a dropdown menu showing "2017_04_11", and "Object:" with a list of four items: "left_femur", "left_tibia", "right_femur", and "right_tibia". Each item has an unchecked checkbox to its left.

- **Subject** This field is automatically populated with the ID of the subject when you load a subject file.
- **Session** This drop down menu is automatically populated with all of the sessions available for the current subject.
- **Object** This list is automatically populated with all of the subject objects when you select a session. If the object has a *Surface File* or *ROI Surface File* specified for it, the object is enabled in the list so you can select it. The object that is selected is displayed in the 3D view. You can select only one object at a time.

From:

<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:

https://has-motion.com/wiki/doku.php?id=other:dsx:orient3d:object_configuration&rev=1721229706

Last update: **2024/07/17 15:21**

