


Surface3D Settings

The Surface3D Settings dialogue allows users to modify some of the application's [parameters](#) as well as the application's color theme.

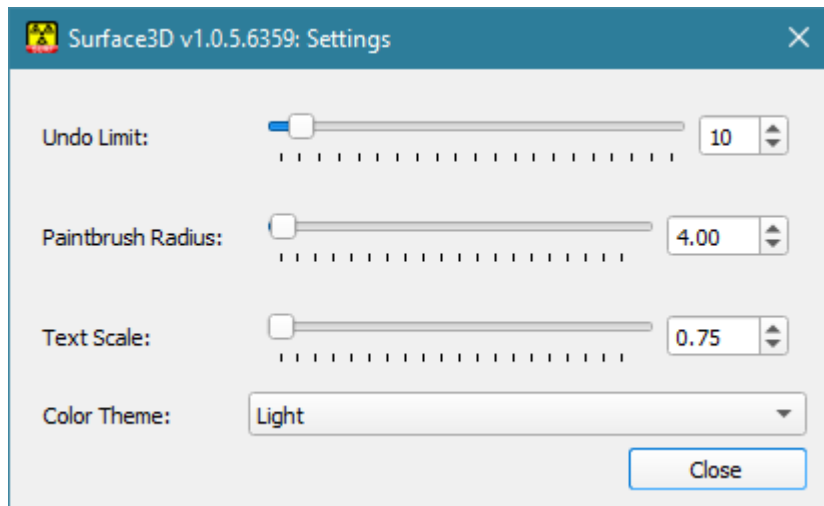

- **Undo Limit** The maximum number of label editing actions that can be undone. The default value is 10.
- **Paintbrush Radius** The radius (in voxels) of the circular paintbrush that is used for painting and erasing segmentation labels. The default value is 5. The radius can also be increased or decreased using the mouse wheel, after pressing **p** to paint or **e** to erase.
- **Text Scale** The scale factor for POI and landmark text labels. The default value is 0.75.
- **Color Theme** You can choose between a light or a dark color theme.

From:
<https://has-motion.com/wiki/> - **Wiki Documentation Site**

Permanent link:
https://has-motion.com/wiki/doku.php?id=other:dsx:surface3d:surface3d_settings

Last update: **2025/05/26 15:44**

