


DSX Configuration

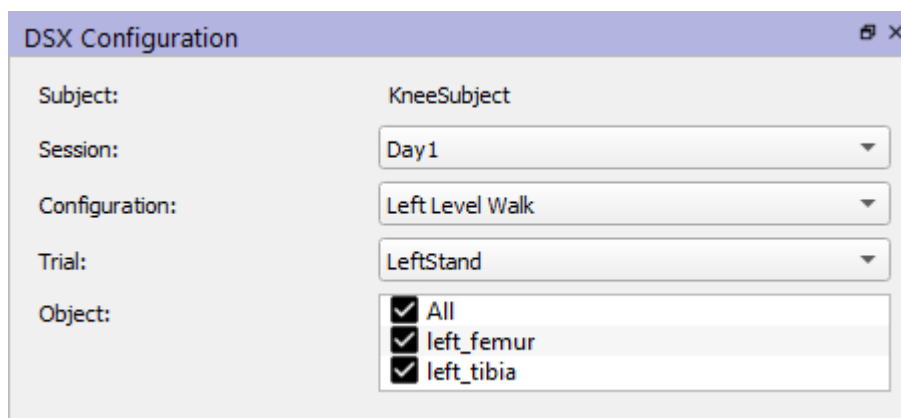

The screenshot shows a window titled "DSX Configuration" with a close button. It contains five fields:

- Subject:** Kneesubject
- Session:** Day1 (dropdown menu)
- Configuration:** Left Level Walk (dropdown menu)
- Trial:** LeftStand (dropdown menu)
- Object:** A list with three items, each with a checked checkbox: All, left_femur, and left_tibia.

- **Subject** This field is automatically populated with the ID of the subject when you load a subject file.
- **Session** This drop down menu is automatically populated with all of the sessions available for the current subject.
- **Configuration** This drop down menu is automatically populated with all of the configurations available for the currently selected session.
- **Trial** This drop down menu is automatically populated with all of the trials available for the currently selected configuration.
- **Object** This list is automatically populated with all of the objects in the currently selected trial. Each object is enabled if it is trackable, meaning that it has a 3D image file associated with it, for generating DRRs.

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
https://has-motion.com/wiki/doku.php?id=other:dsx:x4d:dsx_configuration&rev=1721231062

Last update: 2024/07/17 15:44

