


Settings

X4D's Settings widget allows users to specify how the application appears and behaves.

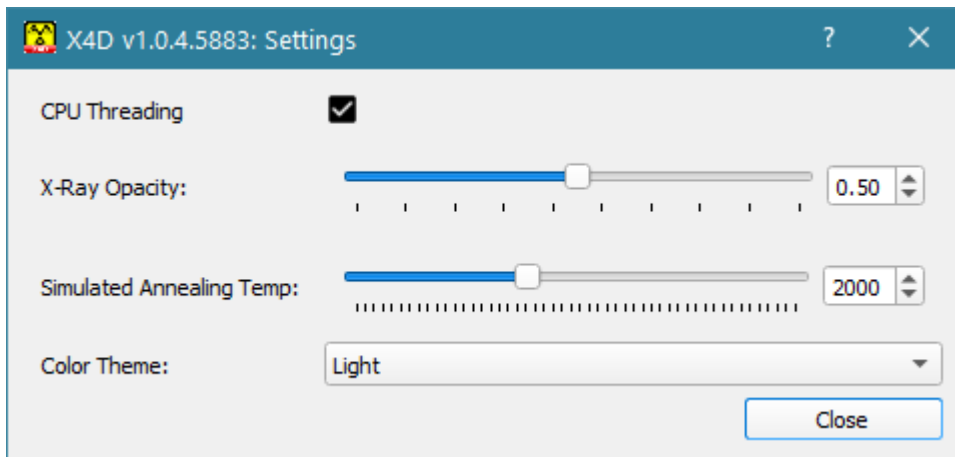

- **CPU Threading** With two CUDA-capable graphics cards, object tracking can be made faster when each card is controlled by its own CPU thread.
- **X-Ray Opacity** The opacity of the X-ray images in the 2D and [3D View](#) windows. Their opacity should be less than 1.0 so that the DRRs can be seen superimposed on them. The default value is 0.5.
- **Simulated Annealing Temp** The initial temperature for the [ASA simulated annealing algorithm](#). The higher the initial temperature, the more slowly the search range will narrow during the optimization. The default value is 2000.
- **Color Theme** You can choose between a light and dark color theme for the graphical interface.

From:
<https://has-motion.com/wiki/> - **Wiki Documentation Site**

Permanent link:
<https://has-motion.com/wiki/doku.php?id=other:dsx:x4d:settings>

Last update: **2025/06/19 19:54**

