

Driveline

Through the [OpenBiomechanics Project](#), [Driveline Baseball](#) has made elite-level motion capture data publicly available. The data can be downloaded from [Driveline's GitHub repository](#).

Having access to large datasets is a great asset, but processing these large datasets can take time (and patience). Inspect3D is designed to streamline this process. The steps in this tutorial describe how to create CMZ files from the pitching and batting data sets, then inspect and group the data in Inspect3D.

Download

Download the files directory from the [repository](#) (Code → Download ZIP).

Data Structure Outline in openbiomechanics-main folder:

- baseball_hitting
 - code
 - v3d
 - CMO.v3s
 - model
 - hitting_v1_model_hybrid.mdh
- data
 - metadata.csv
 - c3d
- baseball_pitching
 - code
 - v3d
 - CMO.v3s
 - model
 - v6_model_hybrid_lm.mdh
 - data
 - metadata.csv
 - c3d

Download the [HAS-Motion Specific Files](#):

1. Filter Script
2. Event Script

1. Note: The event script updates the force platform parameters. Visual3D automatically combines the forces when one foot is in contact with two plates but requires the surfaces to be parallel (or close to). This script adjusts the corners which are coincident to allow the forces to be combined. This script is specific to this data set.

Additional files are provided as examples, but not necessary for processing (files to automatically build group definitions, a subset of group definitions, and sample workspaces).

Outline

1. Build CMZs
 1. This step is done once, or any time you collect more data
2. Build Group definition
3. Inspect data

Build CMZ Files

Pitching

Build CMZ Files for Pitching Data:

1. Click [Load Library](#) Toolbar option
2. In Library Path Dialog, click **Build CMZs**
3. In the Build CMZ Dialog:
 - Set the path to: C:\...\openbiomechanics-main\baseball_pitching\data\c3d
 - [Metadata File](#): * C:\...\openbiomechanics-main\baseball_pitching\data\metadata.csv
 - MDH File: C:\...\openbiomechanics-main\baseball_pitching\code\v3d\model\v6_model_hybrid_lm.mdh
 - Add the following scripts and ensure they are in the correct order:
 - C:\...\CMotion Files\FILTER.v3s
 - C:\...\CMotion Files\Events.v3s
 - C:\...\openbiomechanics-main\baseball_pitching\code\v3d\CMO.v3s
1. Click **Create CMZs**
2. [Metadata](#): Check "Subject Specific" next to playing_level
3. Load Library: Once the CMZ files have been created, in the Library Path Dialog make sure the path is set correctly and click **Load**.

Inspect3D

×

CMZ Directory:

Browse

Generic Builder |

☒ Metadata File:

Browse

☒ MDH File:

Browse

☒

Add Script

Remove Script

Filter.v3s

Events.v3s

CMO.v3s

^

v

Create CMZs

Inspect3D

?

×

Height Units:

Weight Units:

	Header	Description	Subject Specific
1	user	Metric	☒
2	session	Metric	☒
3	session_pitch	Text	☒
4	session_mass_kg	Weight	☒
5	session_height_m	Height	☒
6	age_yrs	Metric	☒
7	playing_level	Text	☒
8	pitch_speed_mph	Metric	☐
9	filename_new	Identifier_Dynamic	☐
10	modelname_new	Identifier_Static	☒


Cancel

Create


Hitting

Driveline Hitting Biomechanics

Build CMZ files for Hitting data:				
-----------------------------------	--	--	--	--

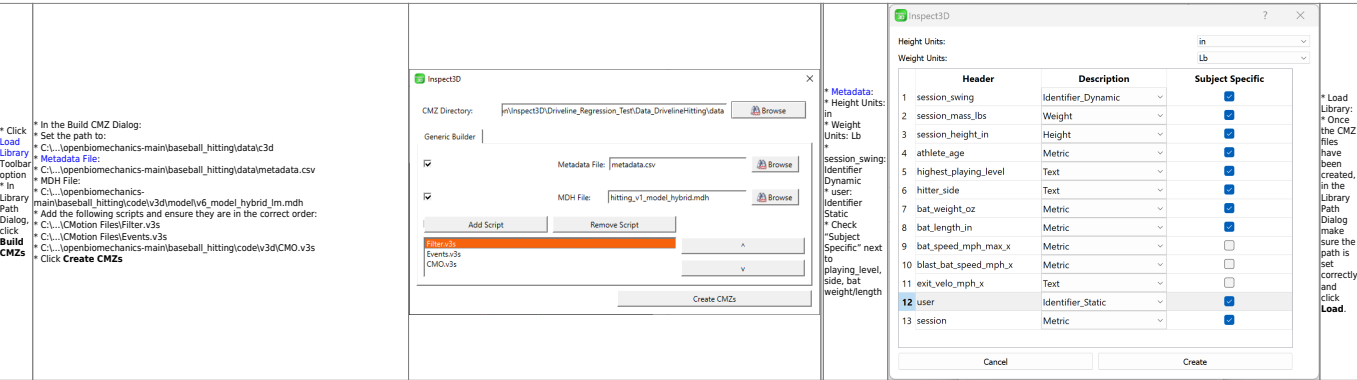

Group Data

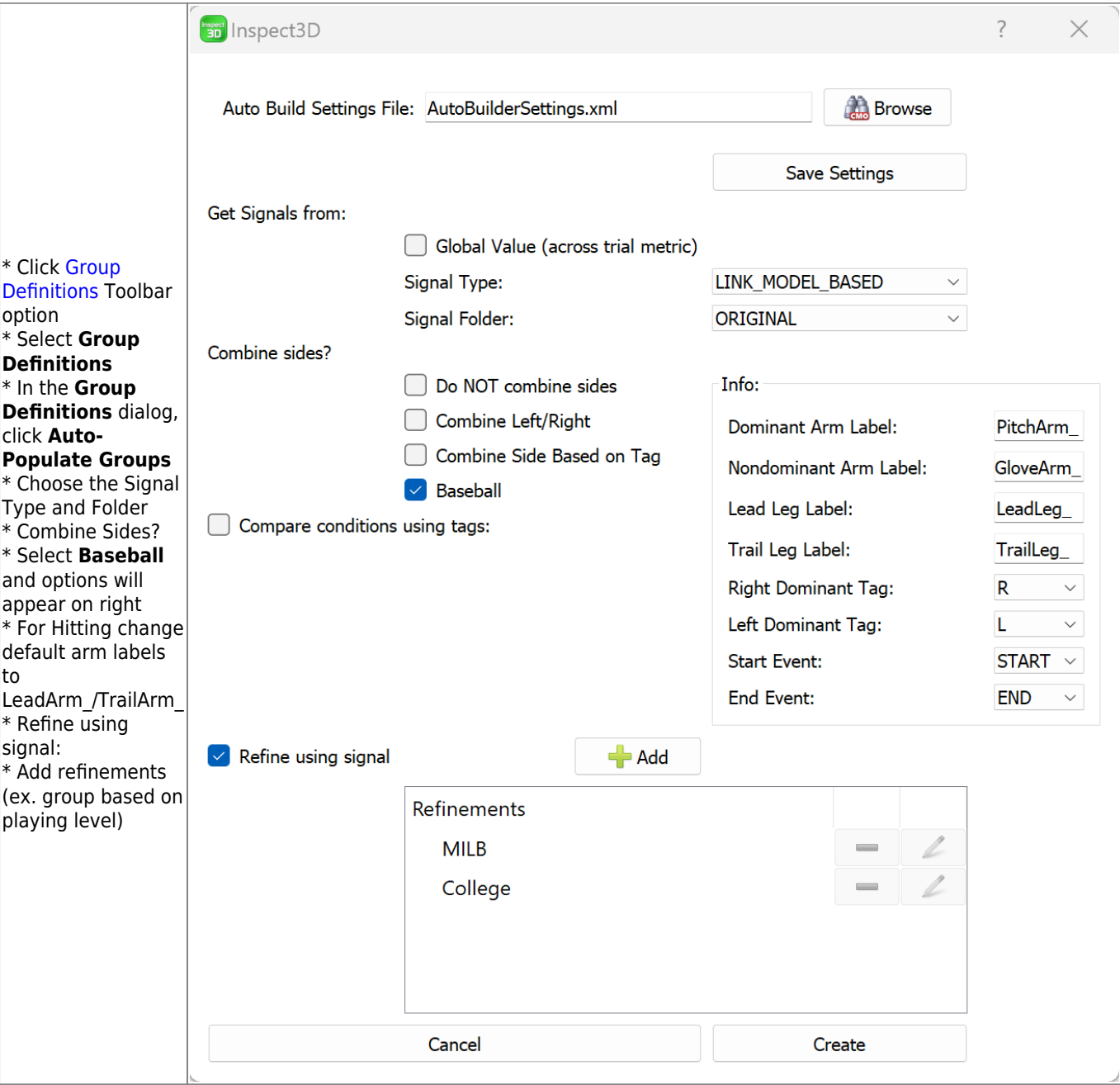

- Alternative: Click **Browse** to load these settings from an XML file to compare two playing levels (MILB vs. College).

- C:\...\CMotion Files\AutoBuilderSettings_Pitch.xml
- C:\...\CMotion Files\AutoBuilderSettings_Hit.xml
- Click **Create**
- In **Group Definition** dialog, delete groups you are not interested in calculating.
- Click **Calculate All Groups**
- Inspect your data!

From:

<https://has-motion.com/wiki/> - **Wiki Documentation Site**

Permanent link:

<https://has-motion.com/wiki/doku.php?id=other:inspect3d:tutorials:driveline&rev=1737057104>

Last update: **2025/01/16 19:51**

