

Excluding Data in Sift

In biomechanics, some trials or trace may not be suitable for analysis. A missed marker, an incorrect force plate assignment, or an unexpected movement during a trial can all introduce noise into the results. Rather than deleting these trials outright, Sift allows users to flag them for exclusion, which keeps the analysis clean while preserving the original data.

Data cleaning is an ongoing process, and exclusions can be revisited at any time as new issues are identified.

How Exclusions Work

Excluding a trace removes it from analysis (and optionally, plots), but does not delete the underlying data. Exclusions can be undone at any time using Ctrl+Z or by right-clicking a workspace and selecting Re-Include Data, and excluded traces can be toggled back on to review them alongside valid data.

Sift supports three levels of exclusion:

- **Individual trace** — exclude a single signal cycle only
- **All traces at those frames** — exclude all queried signals during the same cycle
- **Entire trial** — exclude all data from the affected trial

Writing Exclusions Back to the Data

Once exclusions have been identified, they can be written back to the original workspace files. This is done by adding a “Bad” event to the excluded cycle.... See [excluded_traces](#)

It is also possible to tag entire trials that contain excluded data, making it easy to identify and filter problematic files downstream.

Plot Options

Two display options in the plot are available:

- **Hide Excluded Data:** removes excluded traces from the plot entirely.
- **Show Excluded Data:** keeps excluded traces visible in the plot while still omitting them from analysis.

Tutorials


- [Clean your Data](#) — A step-by-step walkthrough of using Sift to identify and exclude bad trials

across a full dataset, using a ground reaction force quality check as an example.


Example

Excluding Traces Directly


Select the trace(s) to exclude in the Explore Page plot.

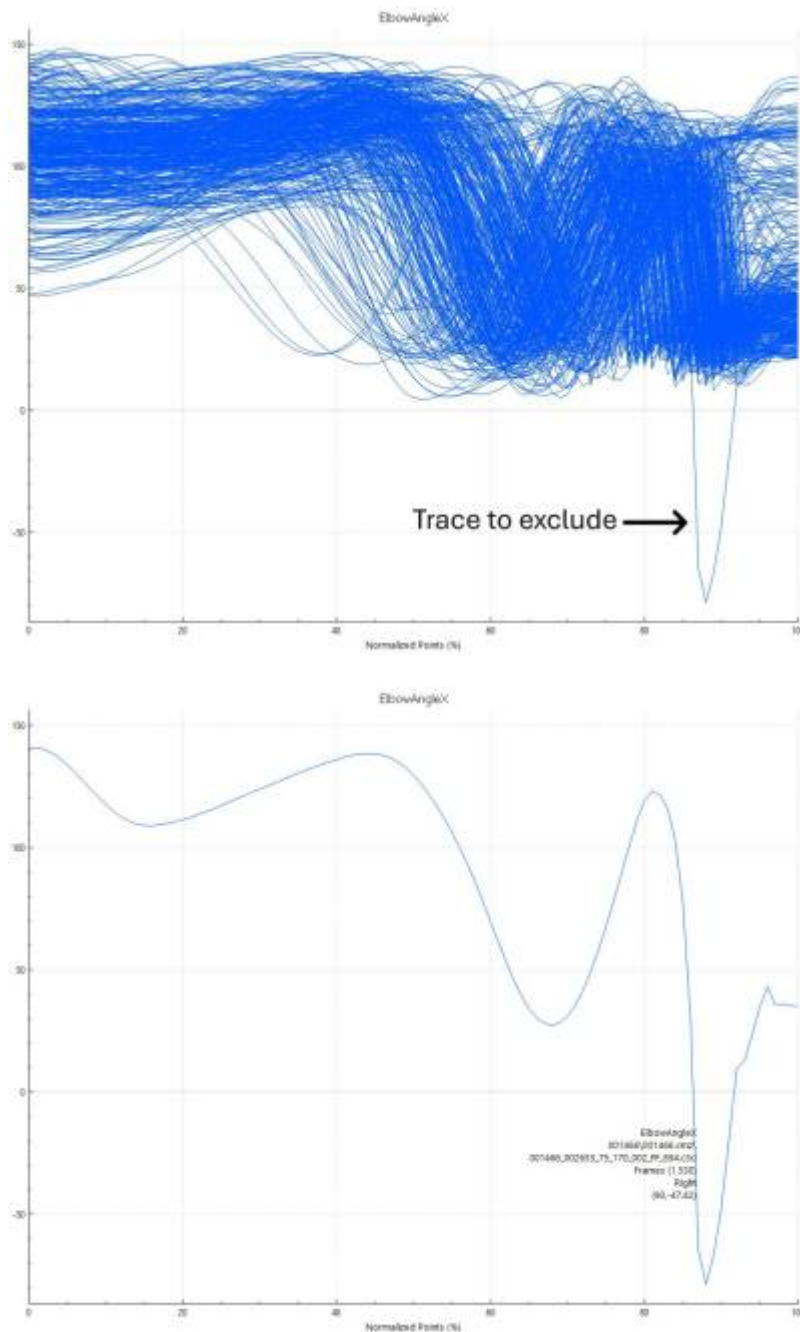

Exclude the selected traces by right-clicking the screen, and selecting **Exclude** → **Exclude Trace (raw data)**, or press **Ctrl+X**

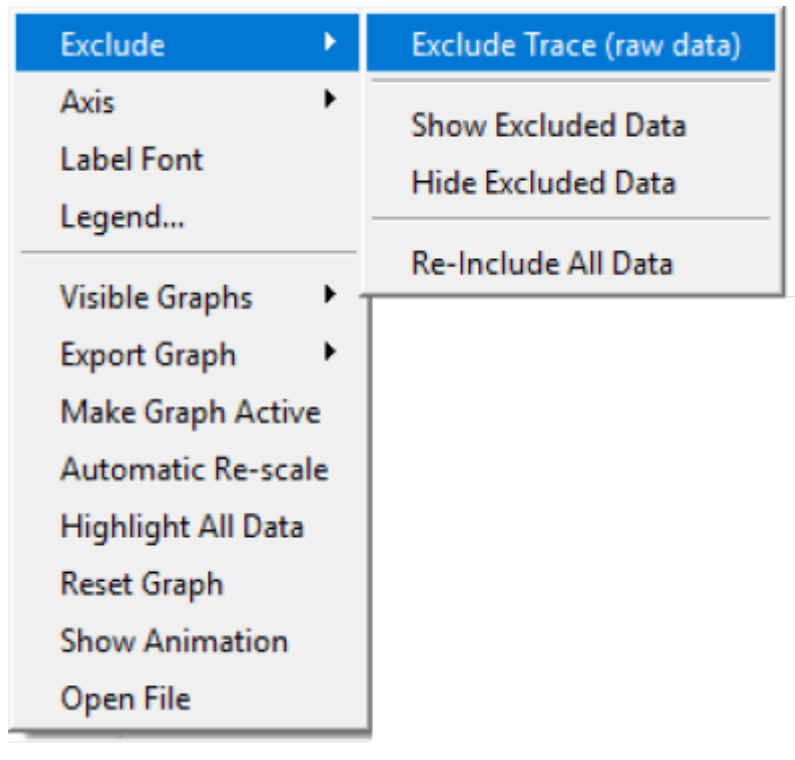

The plot will refresh without the excluded data.

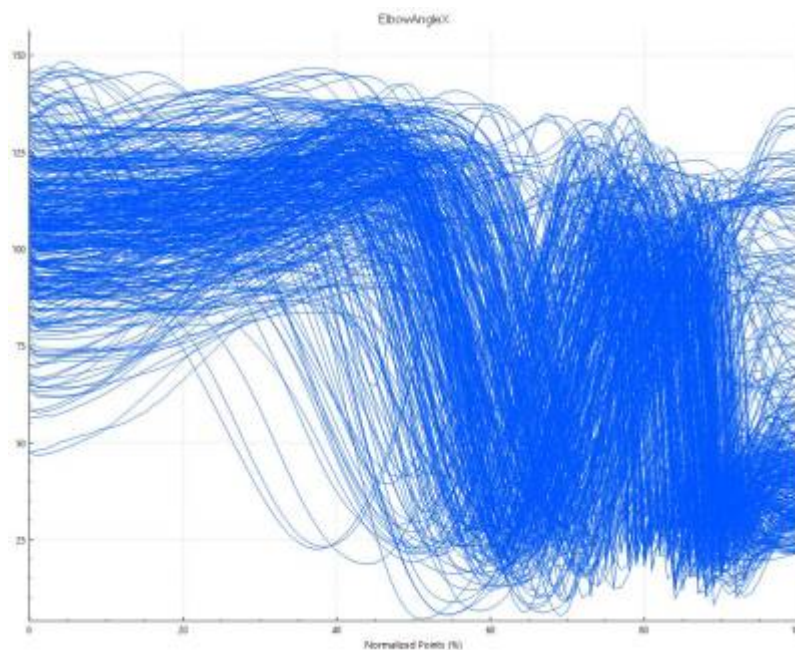

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
https://has-motion.com/wiki/doku.php?id=sift:application:exclude_data&rev=1773154681

Last update: 2026/03/10 14:58

