

Troubleshooting

Listed below are the solutions to common problems that have been reported to the Has-Motion support team. For step-by-step instructions on how to use the various features available in Sift visit our [Tutorials Page](#). If you have a problem that has not been covered below or in the Sift tutorials feel free to e-mail the support team at support@has-motion.ca.

Working Data

All Grouped Data

All data/workspaces that were queried, regardless of what is plotted.

All Plotted Groups

Only perform work on the groups that are plotted on all plots. Work can be performed on:

- All workspaces for those groups (regardless on if they are plotted or not)
- Plotted Workspaces (works on entire workspace, even if only one trace/C3D is plotted)

All Workspaces

Only perform work on the groups that are plotted in the active plot.


- All workspaces for those groups (regardless on if they are plotted or not)
- Plotted Workspaces (only works on the plotted data, so if only one trace/C3D is plotted, that is what will be used)

All Highlighted Data

Only data that is highlighted

Graph is missing data


If you plot data and there are empty data points, you may end up with a straight line at -9999999:

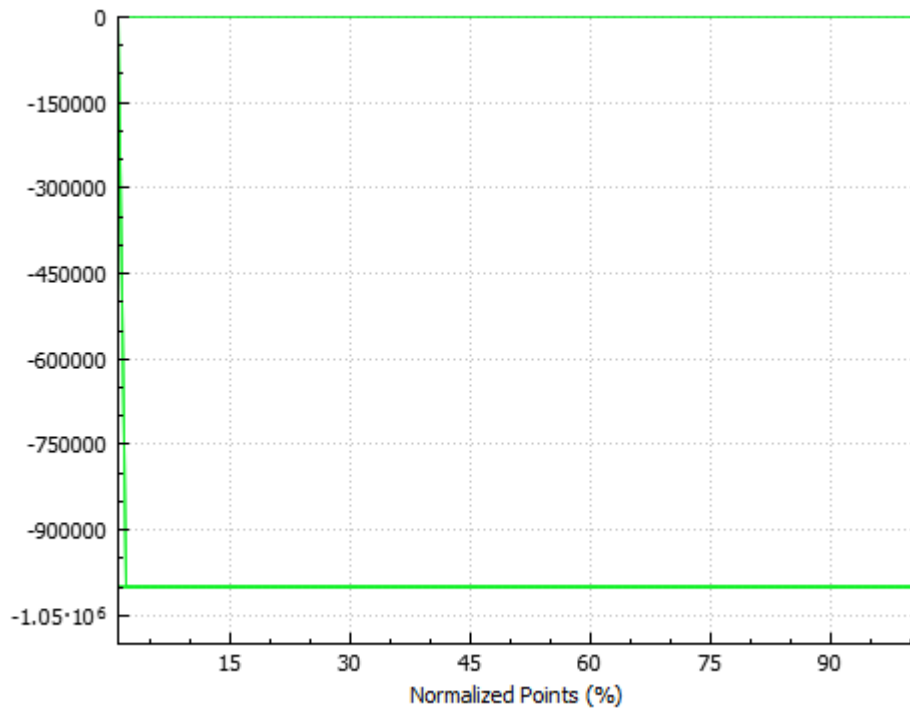

Most calculations (mean, standard deviation, etc) will handle empty data points without an issue, so you can choose not to plot data not found:

- In the [Options Dialogue](#)
- Uncheck **"Show data not found"**
- Click **"Refresh Plot"**

You can also choose to exclude this data by selecting this data, right clicking → Exclude → Exclude trace (raw data).

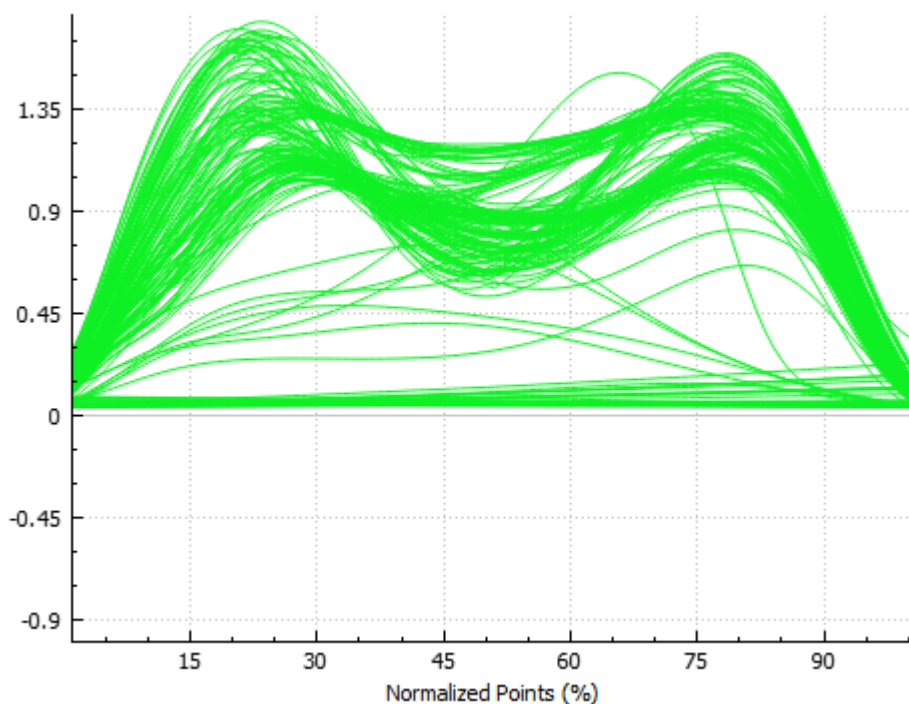

NOTE: Most calculations handle data not found, however, the input to PCA should not contain any missing data points.

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
<https://has-motion.com/wiki/doku.php?id=sift:documentation:troubleshooting&rev=1721230918>

Last update: **2024/07/17 15:41**

