



Create SPM Dialog

The create SPM Dialog is where you complete the process of a [SPM Analysis](#) in Sift. The dialog can be found on the SPM→Mapping Page, and is available after you have created and selected a [GLM](#).

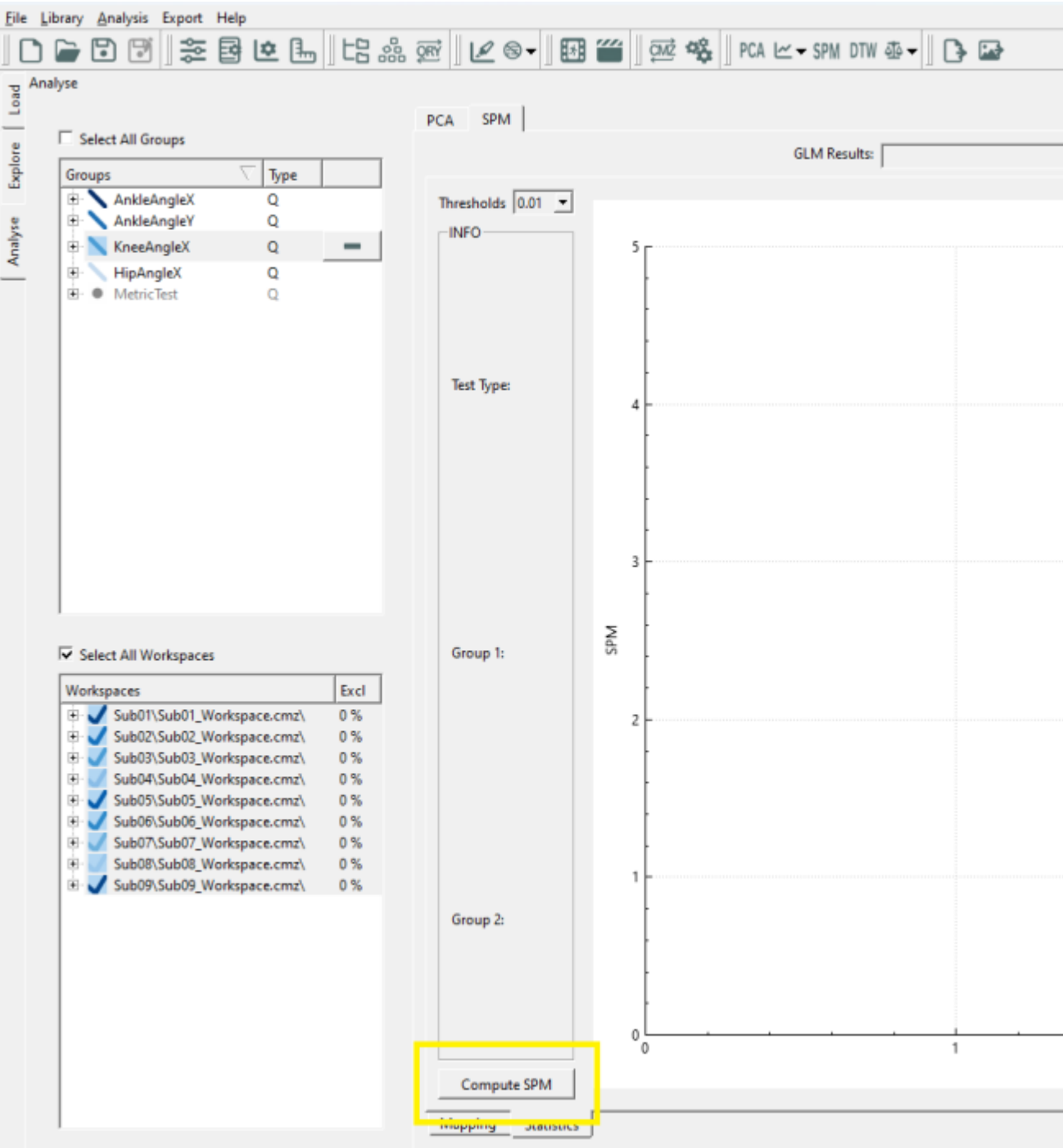

Create SPM Overview

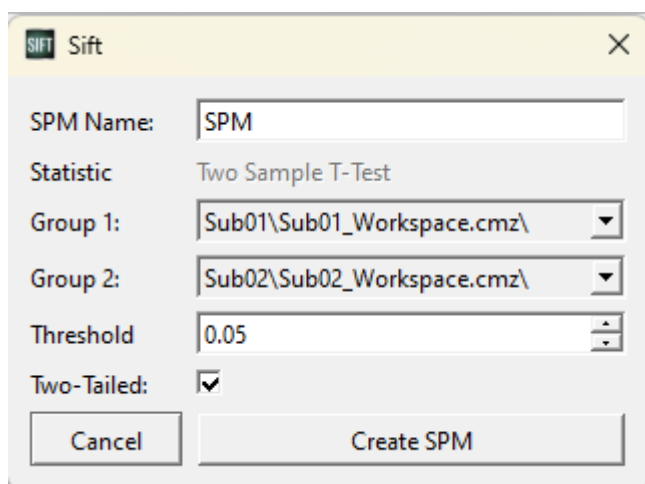

Before open the dialog, ensure you have created and selected a GLM which you would like to create an SPM on. There is currently 1 statistical test to create an SPM with:

- **Two-Sample T-Test:** The standard students two-sample t-test.

Dialog Sections

The dialog has 6 sections:

- **SPM Name:** Enter the name for your SPM . Allows you to cycle through and select different SPM's. Must be a distinct name within the GLM.
- **Statistic:** Select the statistic to create a SPM for.
- **Group 1:** Select the grouping for the first comparison for the t-test. Auto-populated with the groupings for the GLM, whether it is by group or workspace.
- **Group 2:** Select the grouping for the second comparison for the t-test. Auto-populated with the groupings for the GLM, whether it is by group or workspace.
- **Cancel:** Cancels the process of creating a SPM.
- **Create SPM:** Creates a SPM. Makes this the “active” SPM, and populates the [Visualization](#)

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
https://has-motion.com/wiki/doku.php?id=sift:statistical_parametric_mapping:create_spm_dialog&rev=1724857908

Last update: 2024/08/28 15:11

