



Create SPM Dialog

The create SPM Dialog is where you complete the process of a [SPM Analysis](#) in Sift. The dialog can be found on the SPM→Mapping Page, and is available after you have created and selected a [GLM](#).

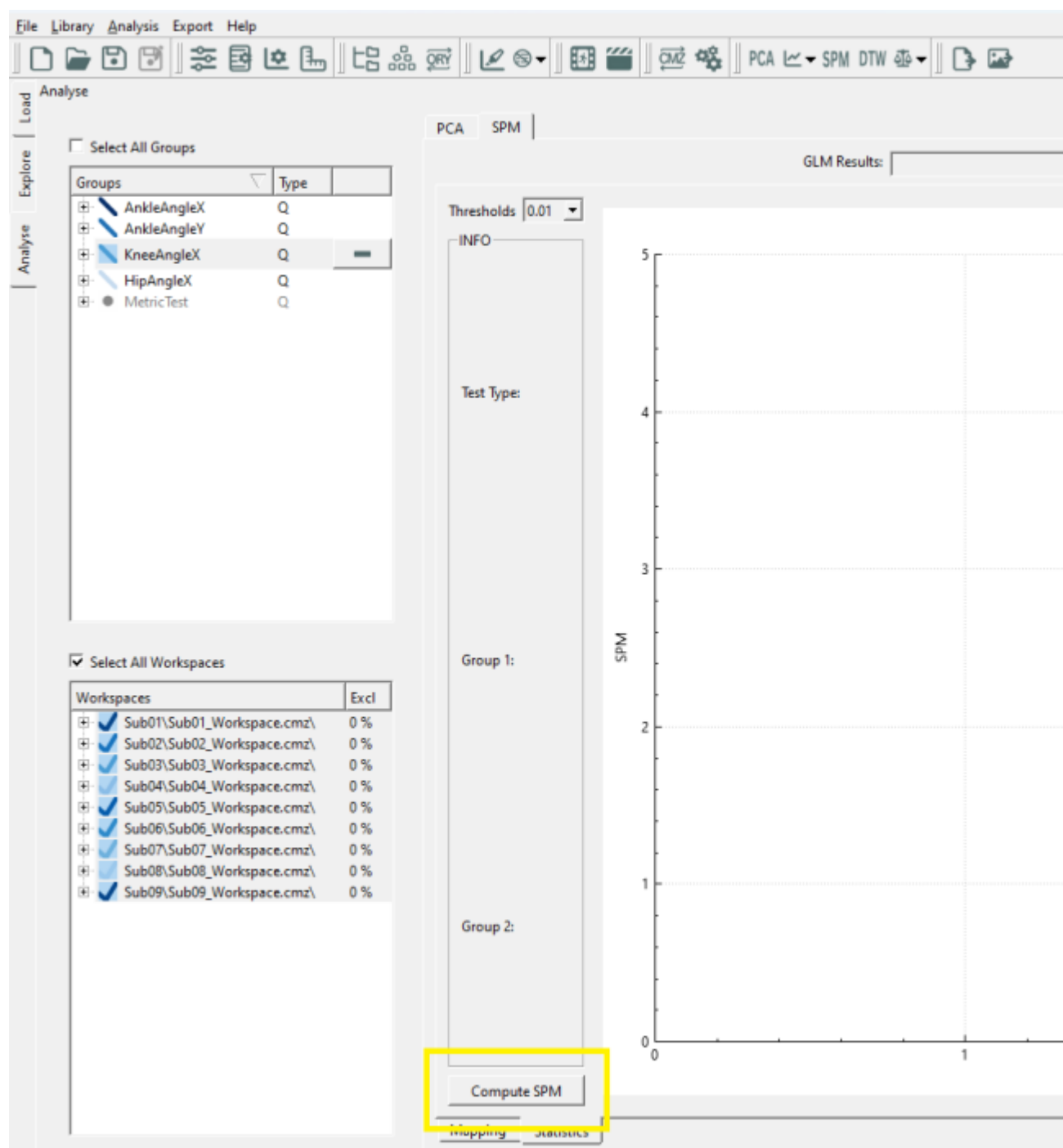

Create SPM Overview

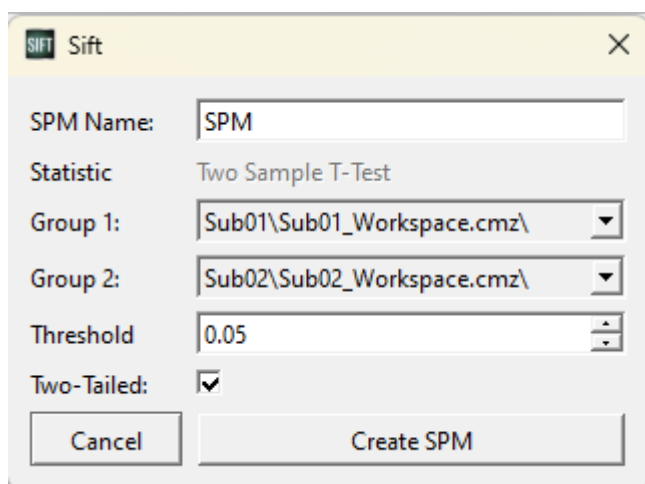

Before open the dialog, ensure you have created and selected a GLM which you would like to create an SPM on. There are currently 4 statistical tests to create an SPM with:

- **One-Sample T-Test:** The standard students one-sample t-test. Tests that the value is non-zero.
- **Two-Sample T-Test:** The standard students two-sample t-test. Tests that the two samples are equivalent.
- **Paired T-Test:** The standard students paired t-test. Tests that all pairs of data are equivalent. Equivalent to a one-sample t-test of the paired differences.
- **ANOVA:** The standard ANOVA test. Tests that multiple groups are equivalent.

Dialog Sections

The dialog has 6 sections:

- **SPM Name:** Enter the name for your SPM . Allows you to cycle through and select different SPM's. Must be a distinct name within the GLM.
- **Statistic:** The statistic that we are testing for. Locked to that of the GLM.
- **Group 1**(If relevant): Select the grouping for the first comparison (ex. for one-sample t-test). Auto-populated with the groupings for the GLM, whether it is by group or workspace.
- **Group 2**(If relevant): Select the grouping for the second comparison (ex. for two-sample t-test). Auto-populated with the groupings for the GLM, whether it is by group or workspace.
- **Threshold **:** **The alpha value to threshold the test against.** * ****Two-Tailed:** If the inference should be two-tailed.
- **Cancel:** Cancels the process of creating a SPM.
- **Create SPM:** Creates a SPM. Makes this the “active” SPM, and populates the [Visualization](#).

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
https://has-motion.com/wiki/doku.php?id=sift:statistical_parametric_mapping:create_spm_dialog&rev=1774966267

Last update: 2026/03/31 14:11

