

Set_Report_Options

```
Set_Report_Options
! /PAGE_LAYOUT=LANDSCAPE
! /PAPER_SIZE=-
! /COLLECTOR=-
! /FACILITY_LINE_2=-
! /ADDRESS_LINE_1=-
! /ADDRESS_LINE_2=-
! /ADDRESS_LINE_3=-
! /PHONE=-
! /FAX=-
! /NOTES_LINE_1=-
! /NOTES_LINE_2=-
! /NOTES_LINE_3=-
! /LOGO=-
! /MARGIN_TOP=-
! /MARGIN_BOTTOM=-
! /MARGIN_LEFT=-
! /MARGIN_RIGHT=-
! /PAGE_NUMBER_POSITION=-
! /PAGE_NUMBER_ALIGNMENT=-
! /PAGE_TAG_LINE_2=-
! /PAGE_TAG_LINE_3=-
! /PAGE_TAG_LINE_4=-
;
```

Manually Set Report Options

These parameters can be changed in the pipeline command, or changed manually by right clicking on the report page and selecting **Report Options**.

Page Layout

Page Layout Tab:

/PAGE_LAYOUT= Portrait or Landscape

/PAPER_SIZE= Letter, Legal or A4

Lab Information

The Lab Information can be added in the Laboratory Information tab.

Laboratory Information Tab:

/COLLECTOR=
/FACILITY_LINE_2=
/ADDRESS_LINE_1=
/ADDRESS_LINE_2=
/ADDRESS_LINE_3=
/PHONE=
/FAX=
/NOTES_LINE_1=
/NOTES_LINE_2=
/NOTES_LINE_3=
/LOGO=


Note that to see the **lab information** in the report, the lab information must be “graphed”.
By default facility line 1 is set to the laboratory name and this parameter cannot be modified.

Page Layout

Page Layout Tab:

1. Margins:
/MARGIN_TOP=
/MARGIN_BOTTOM=
/MARGIN_LEFT=
/MARGIN_RIGHT=
2. Page Number Format:
/PAGE_NUMBER_POSITION=
/PAGE_NUMBER_ALIGNMENT=

Note paper size and orientation can be set in the first two parameters of this pipeline command.

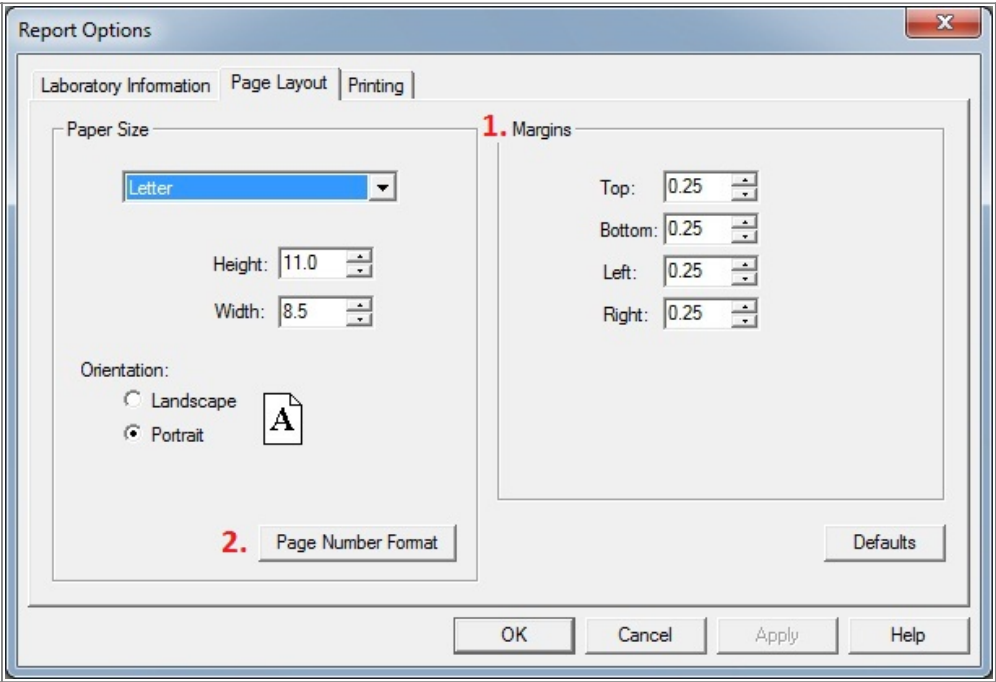

Printing

Printing Tab:

/PAGE_TAG_LINE_2=
/PAGE_TAG_LINE_3=
/PAGE_TAG_LINE_4=

Report Options

Laboratory InformationPage LayoutPrinting

Print

☒Page Titles

☒Page Numbers

☐3-Ring Binder Pips

☐Reverse Print Order

Page Tag Line:

C-Motion TestingPAGE_TAG_LINE_2PAGE_TAG_LINE_3PAGE_TAG_LINE_4

OK

Cancel

Apply

Help

From:
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:
https://has-motion.com/wiki/doku.php?id=visual3d:documentation:pipeline:report_commands:set_report_options&rev=1721229086

Last update: 2024/07/17 15:11

