


Report Mode

Report Dialog

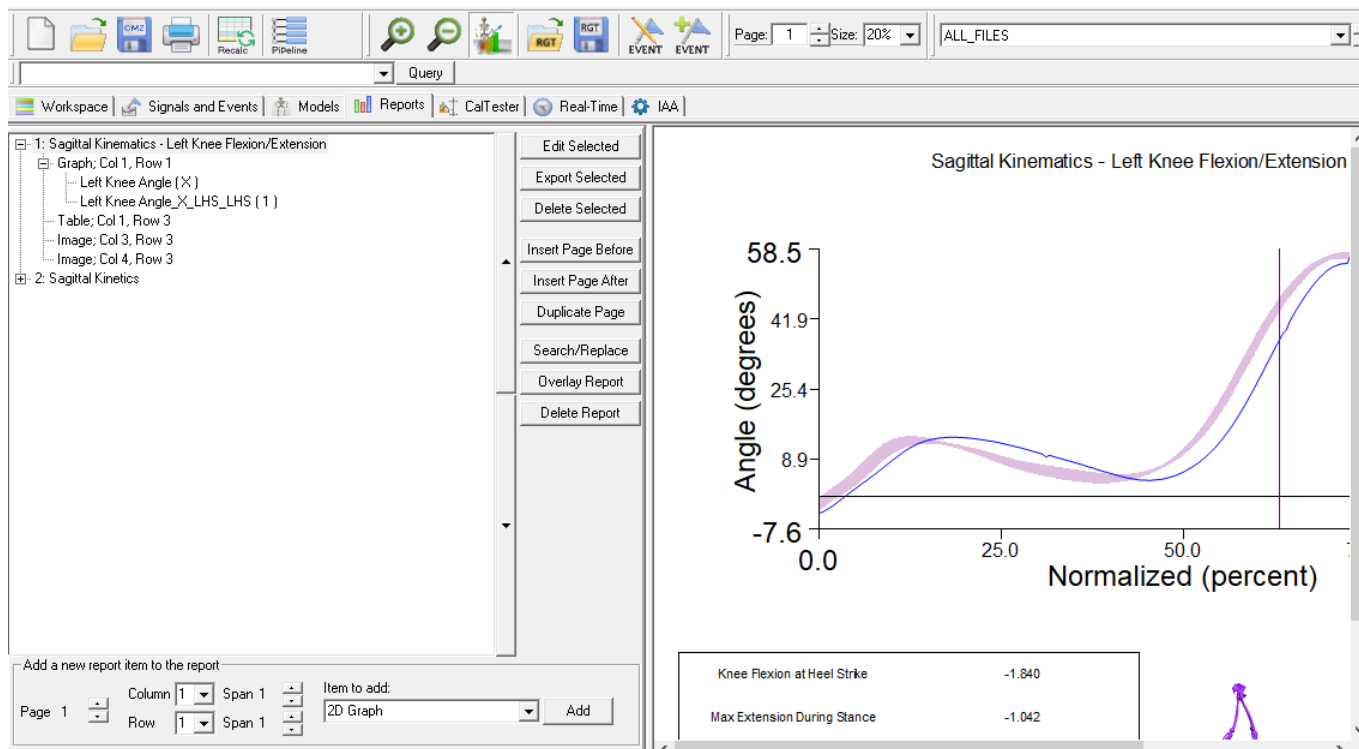

In Visual3Dv6.0#.#, the interface on the Reports tab was modified. The old method of listing the report items in the list was modified to display the report items as a tree control.

The tree control displays the report items as:

- Page Number: Page Title
 - Report Item Type (Graph, Table, Image, Balloon Text, Temporal Distance, Subject Info); Item Location (if the report item has a title (ex. Graph Title is “Knee Flex/Ext”, the graph title will be displayed instead of report item type)
 - Graph Signals (if applicable)


**** ****

Adding a New Item to the Report

New items can be added to the report using the Add New section at the bottom of the dialog.

To add a new graph:

- * If you select a graph from the tree list, and click “Add”, a new graph line will be added to the selected graph
- * If you select a graph from the tree list, and change the graph location (page, column, or row), and click “Add”, a copy of the selected graph will be created at the new location.

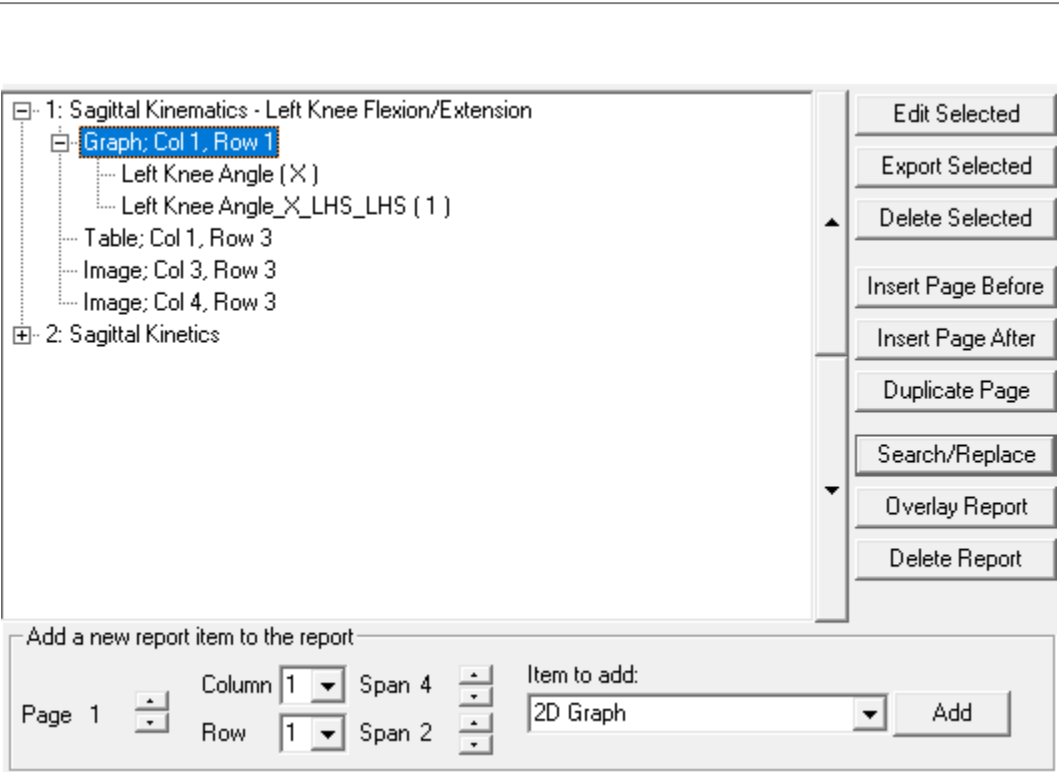


Editing the Report

You can right click on any report item (including report page) to display the context menu:

Right Click gives you the following options:

- * Edit
- * Move
- * Copy (not implemented)
- * Export
- * Delete Selected
- * Page
- * Insert Before
- * Insert After
- * Duplicate
- * Entire Report
- * Search/Replace
- * Overlay
- * Delete Entire Report
- * Add New (not implemented)
- * Modify View

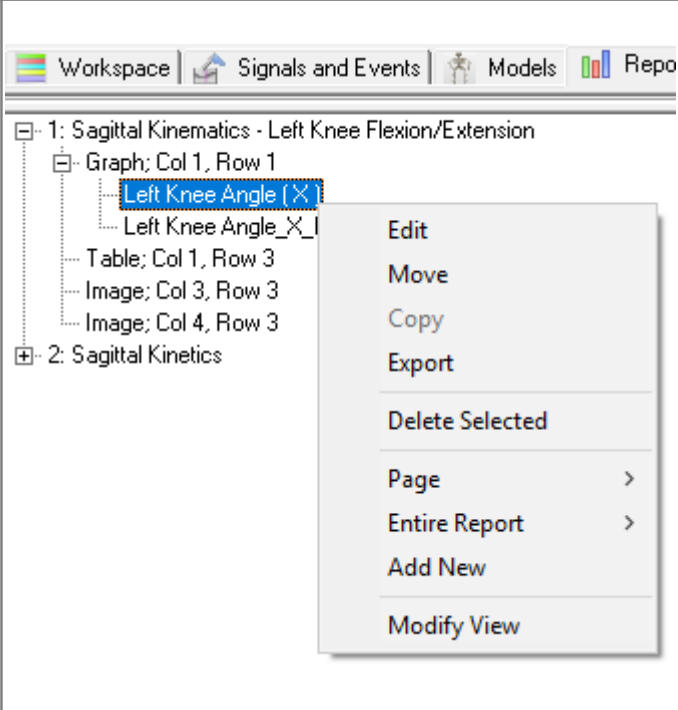

Edit

Brings up the edit dialog for the selected item (page title, graph, etc.)

Move

If you right click on a graph and select move, the graph dialog will appear and the location of the entire graph can be changed. If you right click on a single signal on a graph, a move graph dialog will appear the location of that signal can be changed.

Export

You can export a graph or temporal distance measures from the Report tab using the right click option or by selecting an item and clicking "Export Selected".

Delete Selected

You can delete the selected item (page, graph, table, etc.) by right clicking on an item or selecting the item and clicking "Delete Selected".

Search/Replace

Search/Replace lets you search for a specific item (event label, file tag, or line color) and replace it with another value in the report template.

Overlay

Overlay allows you to quickly overlay a report template - for example if you are plotting a tag "Condition1" you can overlay "Condition2" on top of all graphs plotting "Condition1".

Delete Entire Report

Deletes the entire report - THERE IS NO UNDO BUTTON!

Modify View

Modify View allows you to adjust how the information is displayed in the Report Tab:

- * Show Column Header: will specify what column (1-5) the items are created in and the items in each column will be listed beneath
- * This can be useful if your graphs don't have titles
- * Display Report Edit Buttons: will show/hide report edit buttons
- * This can be useful if you're reviewing a report on a small monitor and not necessarily editing a report
- * For each signal line display (for graphs):
- * Signal name, tag, component, events, etc.
- * This changes what information is displayed in the tree. This can be useful to distinguish between different signals being plotted and for quality assurance to ensure the correct data is being graphed.

Adjust Display

☐ Show Column Header

☒ Display Report Edit Buttons

For each signal line display:

☒ Signal Name

☐ Tag Name

☒ Component

☐ Start Event

☐ End Event

☐ Exclude Events

OK

Cancel

From: <https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link: https://has-motion.com/wiki/doku.php?id=visual3d:documentation:reports:report_mode&rev=1721230967

Last update: 2024/07/17 15:42

