

# Virtual Laboratory

A segment coordinate system in Visual3D does not have to be associated with the subject, or even with motion capture data for that matter. We refer to these as virtual segments. **The only distinction between virtual segments and real segments is that typically we set them to be *Kinematic only* segments, which means that they are not included in the Inverse Dynamics calculations.** In all other ways they are just segments and are constructed with the same rules as all segments.

## Example: Simple Virtual Laboratory

A Virtual Laboratory is not a physical segment in Visual3D but is equivalent to a physical segment to the effect that it has a Segment Coordinate System that can be used to represent the laboratory (or other objects) and has an associated surface model that can be animated as part of the scene. A typical use of the Virtual Laboratory is to represent the floor so that motion capture markers can be projected onto the floor.

**Note:** A typical use of the Virtual laboratory is for creating [Normalized Joint Angles and Normalized Segment Angles](#).

Below is a straightforward virtual lab that has a specified constant orientation.

An example of the Simple Virtual Lab can be downloaded [here](#).

## Create 4 Landmarks

First, create 4 landmarks **Lab\_Origin**, **Lab\_X**, **Lab\_Y**, and **Lab\_Z**. To create these landmarks, go to the **Model** tab and select the **Landmarks** tab.

## 1. Create Lab-Origin:

The screenshot shows the 'Lab-Origin' dialog box within a software interface. The 'Landmarks' tab is selected, and the 'Lab-Origin' sub-tab is active. The 'Landmark Name' field is set to 'Lab-Origin'. Under 'Define Orientation Using:', the 'Starting Point' is set to '(Reference)'. The 'Existing Segment' is set to 'LAB'. Under 'Landmark Offset from Start Point (Reference) or Segment Origin', the 'Offset Using the Following ML/AP/AXIAL Offsets' is selected, with X, Y, and Z values all set to 0.0. The 'Offset by Percent' checkbox is unchecked, and the 'Calibration Only Landmark' checkbox is also unchecked. At the bottom, there are buttons for 'Undo Changes', 'Apply', 'Build Model', and 'Close Tab'.

1. Click **Landmarks** button
2. Click **Add New Landmark** button
3. Create Landmark: *Lab-Origin*

**Landmark Name:** *Lab-Origin*

**Define Orientation Using:**

**Starting Point:** *Leave Blank*

**Existing Segment:** *LAB*

Offset Using the Following **ML/AP/AXIAL** Offsets:

**X:** *0.0* **Y:** *0.0* **Z:** *0.0*

**DO NOT** Check: *Offset by Percent (1.0 = 100%)*

**DO NOT** Check: *Calibration Only Landmark*

Workspace | Signals and Events | Models | Reports | Cal

Segments | Landmarks | Muscles | Subject Data / Metrics

Landmarks | Functional | Digitizing | Lab\_Origin

Landmark Name: Lab\_X

Define Orientation Using:

Starting Point (Reference)

Targets and/or Landmarks:

Ending Point (On a line)

\*Lateral object (On a plane)

\*Project From (Projection onto a line or plane)

\* = Optional

Existing Segment LAB

Landmark Offset from Start Point (Reference) or Segment Origin

Offset to Existing Calibration Target or Landmark

Offset Using the Following ML/AP/AXIAL Offsets

X 0.0

Y 0.0

Z 0.0

Offset by Percent (1.0 = 100%) (Meters when not checked)

Calibration Only Landmark (Not generated for assigned motion file(s))

Undo Changes Apply Build Model Close Tab

## 2. Create Lab\_X:

1. Click **Landmarks** button
2. Click **Add New Landmark** button
3. Create Landmark: *Lab\_X*

**Landmark Name:** *Lab\_X*

**Define Orientation Using:**

**Starting Point:** *Leave Blank*

**Existing Segment:** *LAB*

Offset Using the Following **ML/AP/AXIAL** Offsets:

**X:** *0.1*   **Y:** *0.0*   **Z:** *0.0*

**DO NOT** Check: *Offset by Percent (1.0 = 100%)*

**DO NOT** Check: *Calibration Only Landmark*

Workspace | Signals and Events | Models | Reports | Cal

Segments | Landmarks | Muscles | Subject Data / Metrics

Landmarks | Functional | Digitizing | Lab\_Origin

Landmark Name: Lab\_X

Define Orientation Using:

Starting Point (Reference)

Targets and/or Landmark:

Ending Point (On a line)

\*Lateral object (On a plane)

\*Project From (Projection onto a line or plane)

\* = Optional

Existing Segment LAB

Landmark Offset from Start Point (Reference) or Segment Origin

Offset to Existing Calibration Target or Landmark

Offset Using the Following ML/AP/AXIAL Offsets

X 0.0

Y 0.0

Z 0.0

Offset by Percent (1.0 = 100%) (Meters when not checked)

Calibration Only Landmark (Not generated for assigned motion file(s))

Undo Changes Apply Build Model Close Tab

### 3. Create Lab\_Y:

1. Click **Landmarks** button
2. Click **Add New Landmark** button
3. Create Landmark: *Lab\_Y*

**Landmark Name:** *Lab\_Y*

**Define Orientation Using:**

**Starting Point:** *Leave Blank*

**Existing Segment:** *LAB*

Offset Using the Following **ML/AP/AXIAL** Offsets:

**X:** *0.0* **Y:** *0.1* **Z:** *0.0*

**DO NOT** Check: *Offset by Percent (1.0 = 100%)*

**DO NOT** Check: *Calibration Only Landmark*

Workspace | Signals and Events | Models | Reports | Cal

Segments | Landmarks | Muscles | Subject Data / Metrics

Landmarks | Functional | Digitizing | Lab\_Origin | Lab\_X

Landmark Name: Lab\_Z

Define Orientation Using:

Starting Point (Reference)

Targets and/or Landmarks:

Ending Point (On a line)

\*Lateral object (On a plane)

\*Project From (Projection onto a line or plane)

\* = Optional

Existing Segment LAB

Landmark Offset from Start Point (Reference) or Segment Origin

Offset to Existing Calibration Target or Landmark

Offset Using the Following ML/AP/AXIAL Offsets

X 0.0

Y 0.0

Z 0.0

Offset by Percent (1.0 = 100%) (Meters when not checked)

Calibration Only Landmark (Not generated for assigned motion file(s))

Undo Changes Apply Build Model Close Tab

#### 4. Create Lab\_Z:

1. Click **Landmarks** button
2. Click **Add New Landmark** button
3. Create Landmark: *Lab\_Z*

**Landmark Name:** *Lab\_Z*

**Define Orientation Using:**

**Starting Point:** *Leave Blank*

**Existing Segment:** *LAB*

Offset Using the Following **ML/AP/AXIAL** Offsets:

**X:** *0.0* **Y:** *0.0* **Z:** *0.1*

**DO NOT** Check: *Offset by Percent (1.0 = 100%)*

**DO NOT** Check: *Calibration Only Landmark*

A subset of 3 of these Landmarks can be used to represent the Segment Coordinate System

## Create the Virtual Laboratory Segment

To create the Virtual Laboratory segment, switch to the **Segment** Tab in model builder mode.

1. Create the Virtual Lab Segment:

1. In the Segment Name combo box type **Virtual Lab**

2. Check the **Kinematic Only** Check Box

3. Select the **Create** button

4. In the **Virtual Lab** tab, enter these values:

**Define Proximal Joint and Radius**

Lateral: None

Joint: Lab\_Z

Medial: None

Radius: 0.01

**Define Distal Joint and Radius**

Lateral: None

Joint: Lab\_Origin

Medial: None

Radius: 0.01

**Extra Target to define Orientation**

Select: Lateral Lab\_X

**Select Tracking Targets:**

Check: Use Calibration Targets for Tracking

5. Click on **Build Model**.

6. Click on **Close Tab** before proceeding.

Workspace | Signals and Events | Models | Reports | Cal

Segments | Landmarks | Muscles | Subject Data / Metrics

Segments | Segment Properties | IK Constraints

Segment Name  
Virtual Lab

Segment Type  
Visual 3D

☒ Kinematic Only

Create

Workspace | Signals and Events | Models | Reports | CalTester | Real-Time | IAA

Segments | Landmarks | Muscles | Subject Data / Metrics

Segments | Segment Properties | IK Constraints | Virtual Lab

Define Proximal Joint and Radius

Lateral: None

Joint Center: Lab\_Z

Medial: None

Radius (Meters): 0.01

Define Distal Joint and Radius

Lateral: None

Joint Center: Lab\_Origin

Medial: None

Radius (Meters): 0.01

Extra Target to Define Orientation (if needed)

Location: Lateral

Lab\_X

Select Tracking Targets (Ctrl-Left Mouse Click to Multiselect)

☒ Use Calibration Targets for Tracking

Depth (Meters):

Close Tab

Guess Properties

Apply

Build Model

Visual3D v6 Prof

Z

## Adding a Surface Model

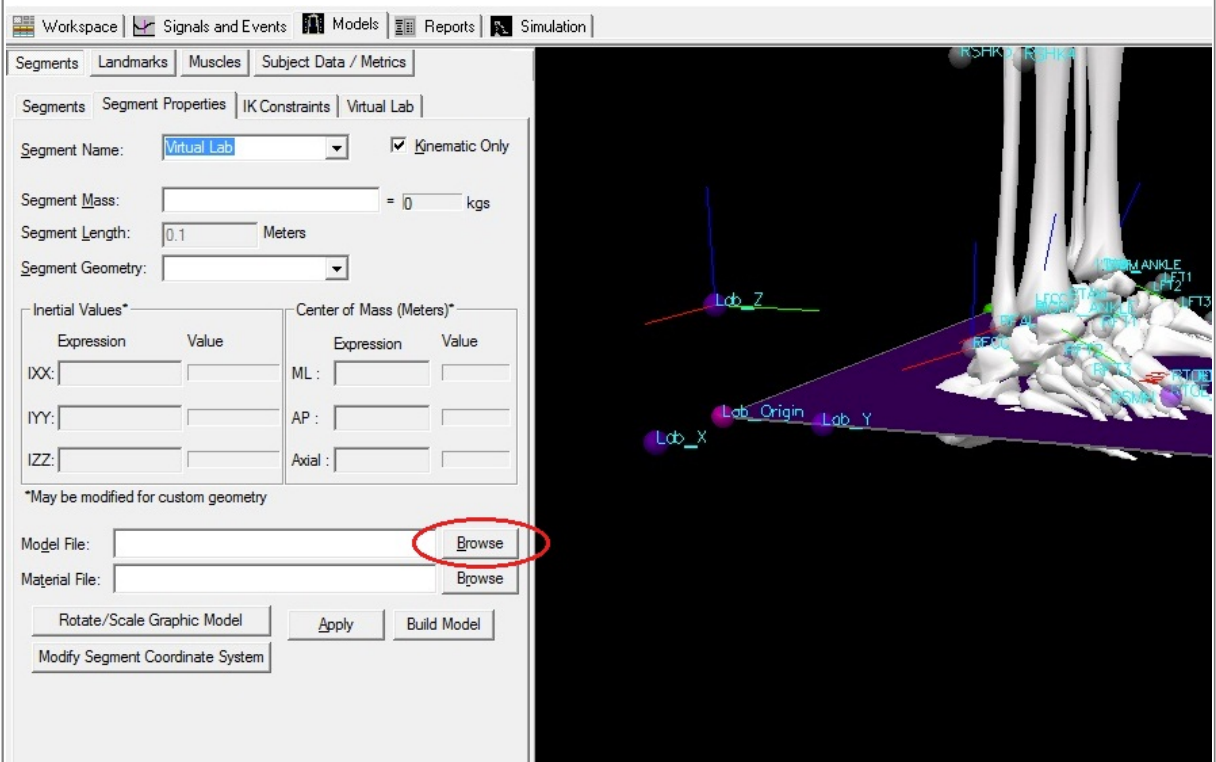
<https://has-motion.com/wiki/>

Printed on 2026/09/15 04:22

To add a surface model or graphic switch to the **Segment Properties** Tab. Select **Virtual Lab** as Segment Name. Near the bottom of the dialog **Browse** for a Model File (.obj or .v3g) of your choosing.

Select **Apply**.

**Note:** A surface model need not be applied to the virtual lab.



## Example: Virtual segment aligned with a different principal axis of the laboratory

The following example will construct a segment coordinate system that is rotated from the laboratory coordinate system.

An example of a Rotated Virtual Lab can be downloaded [here](#).

## Create 3 Landmarks

First, create 3 landmarks **Lab\_Origin**, **Lab\_Other**, and **Lab\_Z**. To create these landmarks, go to the **Model** tab and select the **Landmarks** tab.

|                                                                                                                                                                                        |                                                                                                                                                 |                                                                                                                                                                                                           |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                            |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <div><div>1. Create Lab-Origin:</div><div><div><div>• Click Landmarks button</div><div>• Click Add New Landmark button</div><div>• Create Landmark: Lab-Origin</div></div></div></div> | <div><div>Landmark Name: Lab-Origin</div><div>Define Orientation Using: Starting Point: Leave Blank</div><div>Existing Segment: LAB</div></div> | <div><div>• Offset Using the Following ML/AP/AXIAL Offsets: X: 0.0 Y: 0.0 Z: 0.0</div><div>• DO NOT Check: Offset by Percent (1.0 = 100%)</div><div>• DO NOT Check: Calibration Only Landmark</div></div> | <div><div>Workspace   Signals and Events   Models   Reports   Cal</div><div>Segments   Landmarks   Muscles   Subject Data / Metrics</div><div>Landmarks   Functional   Digitizing   Lab-Origin</div><div>Landmark Name: Lab-Origin</div><div>Define Orientation Using:<div>Starting Point (Reference)</div><div>Targets and/or Landmarks:<div>Ending Point (On a line)</div><div>*Lateral object (On a plane)</div><div>*Project From (Projection onto a line or plane)</div><div>* = Optional</div></div><div>Existing Segment LAB</div></div><div>Landmark Offset from Start Point (Reference) or Segment Origin<div>Offset to Existing Calibration Target or Landmark</div><div>Offset Using the Following ML/AP/AXIAL Offsets<div>X 0.0</div><div>Y 0.0</div><div>Z 0.0</div><div>Offset by Percent (1.0 = 100%) (Meters when not checked)</div></div></div><div>Calibration Only Landmark (Not generated for assigned motion file(s))</div><div>Undo Changes   Apply   Build Model   Close Tab</div></div>            |
| <div><div>2. Create Lab-Other:</div><div><div>• Click Landmarks button</div><div>• Click Add New Landmark button</div><div>• Create Landmark: Lab-Other</div></div></div>              | <div><div>Landmark Name: Lab-Other</div><div>Define Orientation Using: Starting Point: Leave Blank</div><div>Existing Segment: LAB</div></div>  | <div><div>• Offset Using the Following ML/AP/AXIAL Offsets: X: 0.0 Y: 0.1 Z: 0.0</div><div>• DO NOT Check: Offset by Percent (1.0 = 100%)</div><div>• DO NOT Check: Calibration Only Landmark</div></div> | <div><div>Workspace   Signals and Events   Models   Reports   Cal</div><div>Segments   Landmarks   Muscles   Subject Data / Metrics</div><div>Landmarks   Functional   Digitizing   Lab-Origin   Lab-Other</div><div>Landmark Name: Lab-Other</div><div>Define Orientation Using:<div>Starting Point (Reference)</div><div>Targets and/or Landmarks:<div>Ending Point (On a line)</div><div>*Lateral object (On a plane)</div><div>*Project From (Projection onto a line or plane)</div><div>* = Optional</div></div><div>Existing Segment LAB</div></div><div>Landmark Offset from Start Point (Reference) or Segment Origin<div>Offset to Existing Calibration Target or Landmark</div><div>Offset Using the Following ML/AP/AXIAL Offsets<div>X 0.0</div><div>Y 0.1</div><div>Z 0.0</div><div>Offset by Percent (1.0 = 100%) (Meters when not checked)</div></div><div>Calibration Only Landmark (Not generated for assigned motion file(s))</div><div>Undo Changes   Apply   Build Model   Close Tab</div></div></div> |

|                                                                                                                                                                                                     |                                                                                                                                                              |                                                                                                                                                                                                                                                                             |  |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--|
| <p><b>3. Create Lab_Z:</b></p> <ul style="list-style-type: none"> <li>Click <b>Landmarks</b> button</li> <li>Click <b>Add New Landmark</b> button</li> <li>Create Landmark: <b>Lab_Z</b></li> </ul> | <p><b>Landmark Name:</b><br/><b>Lab_Z</b></p> <p><b>Define Orientation Using:</b><br/><b>Starting Point:</b><br/><b>Existing Segment:</b><br/><b>LAB</b></p> | <ul style="list-style-type: none"> <li>Offset Using the Following <b>ML/AP/AXIAL</b> Offsets:<br/><b>X: 0.0 Y: 0.0 Z: 0.1</b></li> <li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li> <li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li> </ul> |  |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--|

## Create the Virtual Laboratory Segment

To create the Virtual Laboratory segment, switch to the **Segment** Tab in model builder mode and enter the following:

|                                                                                                                                                                                                                                                      |  |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--|
| <p><b>1. Create the Virtual Lab Segment:</b></p> <ol style="list-style-type: none"> <li>In the Segment Name combo box type <b>Virtual Lab</b></li> <li>Check the <b>Kinematic Only</b> Check Box</li> <li>Select the <b>Create</b> button</li> </ol> |  |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--|

4. In the **Virtual Lab** tab, enter these values:

**Define Proximal Joint and Radius**

**Lateral:**  
None

**Joint:**  
Lab\_Z

**Medial:**  
None

**Radius:** 0.01

**Define Distal Joint and Radius**

**Lateral:**  
None

**Joint:**  
Lab-Origin

**Medial:**  
None

**Radius:** 0.01

**Extra Target to define Orientation**

**Select:**  
Lateral   Lab\_Other

**Select Tracking Targets:**  
  
**Check:** Use Calibration Targets for Tracking

5. Click on **Build Model**.

6. Click on **Close Tab** before proceeding.

Workspace | Signals and Events | Models | Reports | Simulation

Segments | Landmarks | Muscles | Subject Data / Metrics

Segments | Segment Properties | IK Constraints | Virtual Lab

Define Proximal Joint and Radius

LateralJointMedialRadius (Meters)

NoneLab\_ZNone0.01

Define Distal Joint and Radius

LateralJointMedialRadius (Meters)

NoneLab-OriginNone0.01

Extra Target to Define Orientation (if needed)

LocationLateralLab\_Other

Select Tracking Targets (Ctrl-Left Mouse Click to Multiselect)

☒ Use Calibration Targets for Tracking

Depth (Meters):

Close Tab

Guess Properties

Apply

Build Model

We have just created a virtual lab that is aligned with a different axis of the lab coordinate system (rotated).

## Example: Creating a virtual laboratory that changes with the direction of walking

In this example file (Example3.cmo) the subject is walking in the anterior direction in one walking trial and walking in the posterior direction of the laboratory coordinate system (LAB) in the other trial. If the pelvis angle is defined as the pelvis with respect to the laboratory, the data are rotated 180 degrees from the desired angle for one of the trials (see *Pelvis Angle wrt Lab*). To create a pelvis angle that makes intuitive sense (e.g. in the direction of walking), a virtual lab segment was created. The pelvis angle can then be defined with respect to the virtual lab (see *Pelvis Angle wrt Virtual Lab*). In Example2.cmo the virtual laboratory was defined statically. But, in this example (Example3.cmo), the Virtual lab is defined relative to the lateral direction of the pelvis segment; the lateral direction is forced to lie along one of the principal axes rather than to float with the pelvis.

An example of a Virtual Lab which changes with the direction of walking can be downloaded [here](#).

See Also [A virtual lab with 4 possible orientations](#)

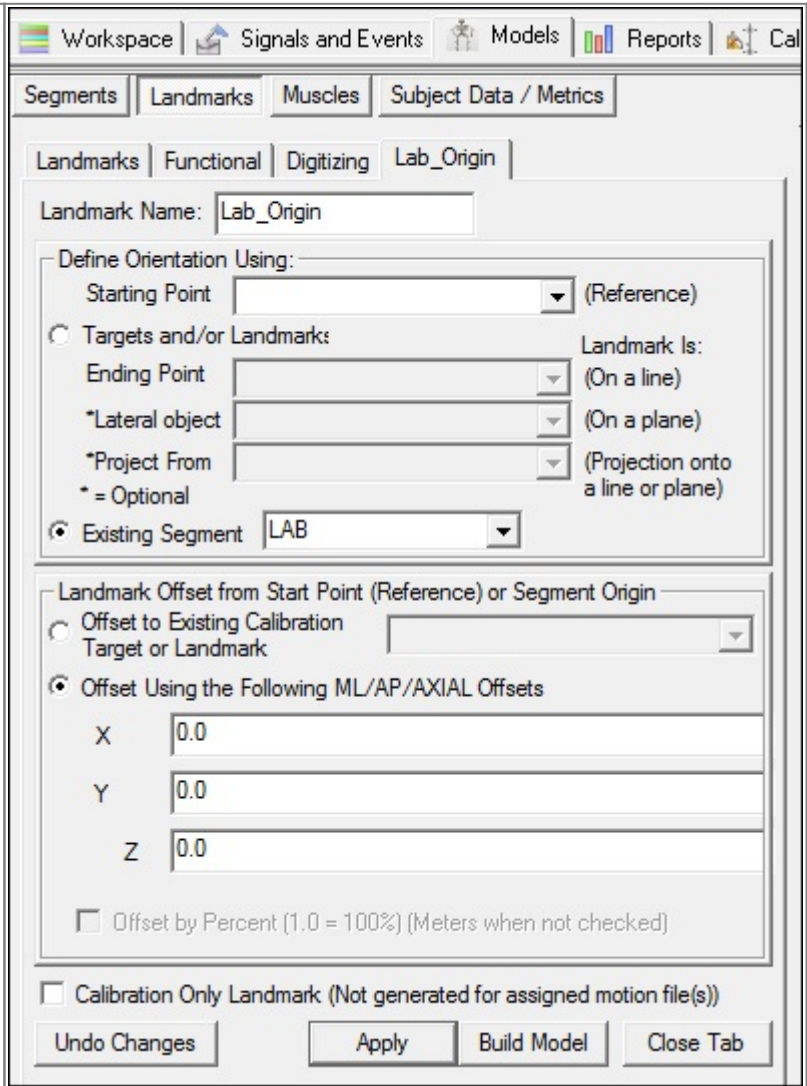
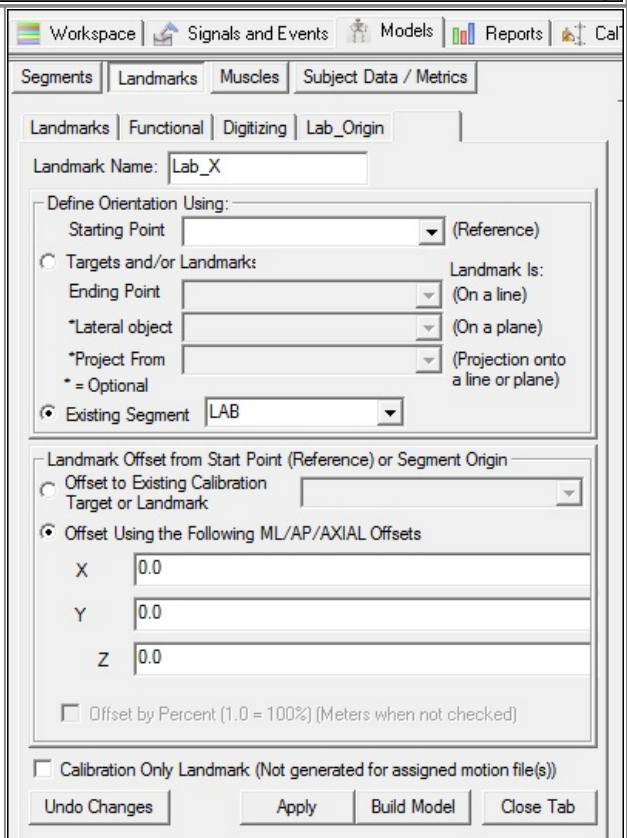
### Create 3 Landmarks

Create 2 landmarks **Lab-Origin**, **Lab-X**, and **Lab-Z**. To create these landmarks, go to the **Model** tab

https://has-motion.com/wiki/

Printed on 2026/09/15 04:22

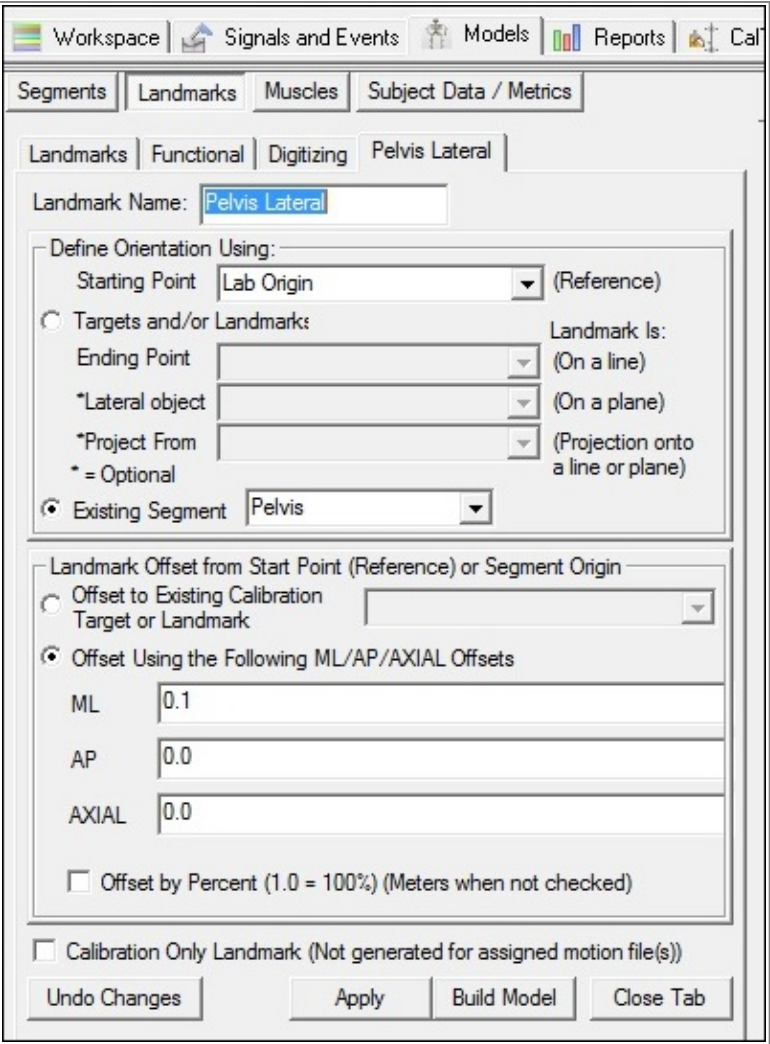
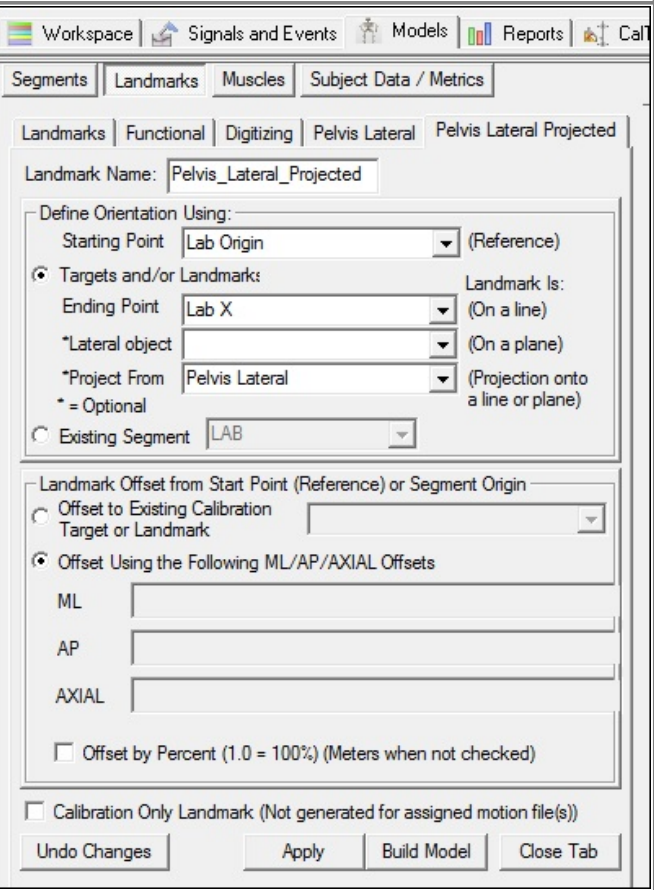
and select the **Landmarks** tab.

|                                                                                                                                                                                                               |                                                                                                                                                                                             |                                                                                                                                                                                                                                                                                           |                                                                                      |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
| <p><b>1. Create Lab_Origin:</b></p> <ul style="list-style-type: none"> <li>Click <b>Landmarks</b> button</li> <li>Click <b>Add New Landmark</b> button</li> <li>Create Landmark: <i>Lab_Origin</i></li> </ul> | <p><b>Landmark Name:</b><br/><i>Lab_Origin</i></p> <p><b>Define Orientation Using:</b><br/><b>Starting Point:</b><br/><i>Leave Blank</i></p> <p><b>Existing Segment:</b><br/><i>LAB</i></p> | <ul style="list-style-type: none"> <li>Offset Using the Following <b>ML/AP/AXIAL</b> Offsets:<br/><b>X:</b> 0.0 <b>Y:</b> 0.0 <b>Z:</b> 0.0</li> <li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li> <li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li> </ul> |   |
| <p><b>2. Create Lab_X:</b></p> <ul style="list-style-type: none"> <li>Click <b>Landmarks</b> button</li> <li>Click <b>Add New Landmark</b> button</li> <li>Create Landmark: <i>Lab_X</i></li> </ul>           | <p><b>Landmark Name:</b><br/><i>Lab_Other</i></p> <p><b>Define Orientation Using:</b><br/><b>Starting Point:</b><br/><i>Leave Blank</i></p> <p><b>Existing Segment:</b> <i>LAB</i></p>      | <ul style="list-style-type: none"> <li>Offset Using the Following <b>ML/AP/AXIAL</b> Offsets:<br/><b>X:</b> 0.1 <b>Y:</b> 0.0 <b>Z:</b> 0.0</li> <li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li> <li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li> </ul> |  |

|                                                                                                                                                                   |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <div><div>3. Create Lab_Z:</div><div><div>• Click Landmarks button</div><div>• Click Add New Landmark button</div><div>• Create Landmark: Lab_Z</div></div></div> | <div><div>Landmark Name: Lab_Z</div><div>Define Orientation Using:</div><div>Starting Point: <div></div> (Reference)</div><div><div>• Targets and/or Landmarks:</div><div>Ending Point: <div></div> (On a line)</div><div>*Lateral object: <div></div> (On a plane)</div><div>*Project From: <div></div> (Projection onto a line or plane)</div><div>* = Optional</div></div><div><div>• Existing Segment: LAB</div></div><div>Landmark Offset from Start Point (Reference) or Segment Origin</div><div><div>• Offset to Existing Calibration Target or Landmark: <div></div></div><div>• Offset Using the Following ML/AP/AXIAL Offsets</div><div>X: <div>0.0</div></div><div>Y: <div>0.0</div></div><div>Z: <div>0.0</div></div><div><input type="checkbox"/> Offset by Percent (1.0 = 100%) (Meters when not checked)</div></div><div><input type="checkbox"/> Calibration Only Landmark (Not generated for assigned motion file(s))</div><div><div>Undo Changes</div><div>Apply</div><div>Build Model</div><div>Close Tab</div></div></div> |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

Create a Landmark from the Pelvis lateral axis projected onto a Laboratory Axis

Create a landmark **Pelvis\_Lateral** 10 cm from the Laboratory Origin in the lateral direction of the Pelvis. Create another landmark **Pelvis\_Lateral\_Projected** that projects **Pelvis\_Lateral** onto a line formed by the **Lab-Origin** and **Lab\_X** landmarks. Then create a landmark

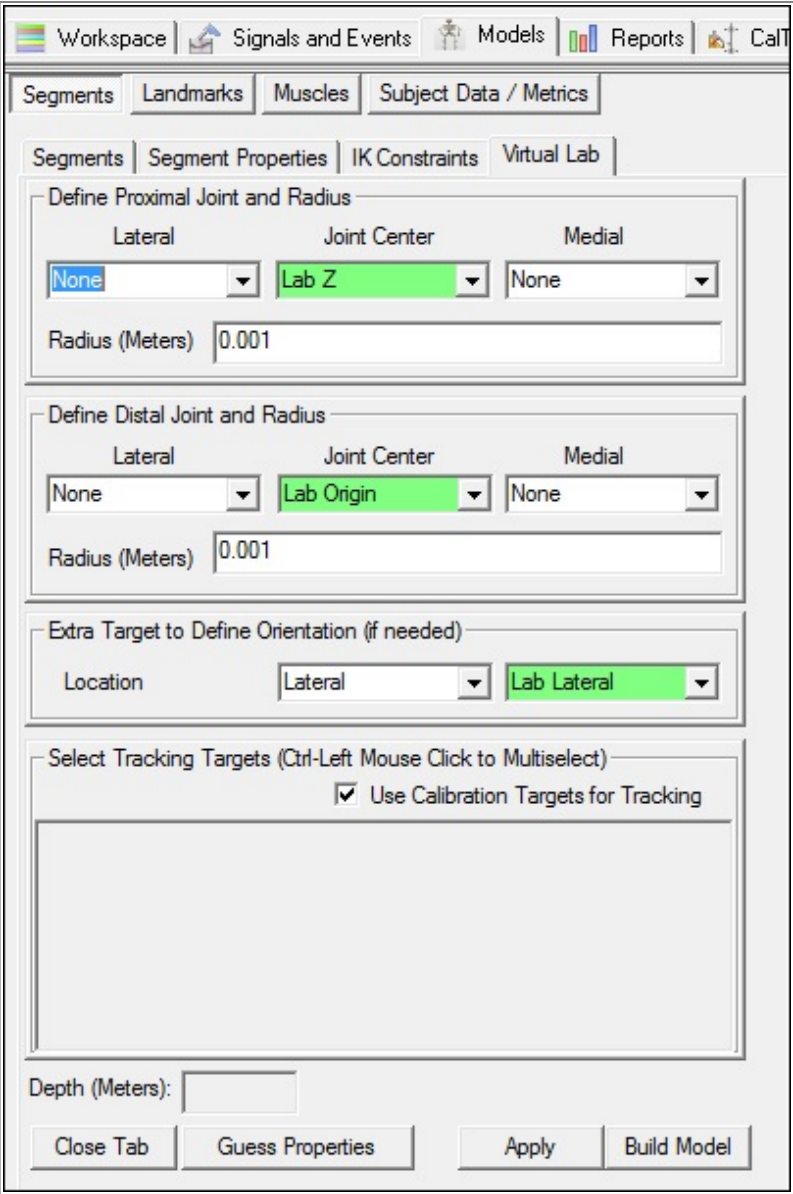
|                                                                                                                                                                                                                                       |                                                                                                                                                                                                                                          |                                                                                                                                                                                                                                                                                               |                                                                                      |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|
| <p><b>1. Create Pelvis_Lateral:</b></p> <ul style="list-style-type: none"> <li>Click <b>Landmarks</b> button</li> <li>Click <b>Add New Landmark</b> button</li> <li>Create Landmark: <i>Pelvis_Lateral</i></li> </ul>                 | <p><b>Landmark Name:</b><br/><i>Pelvis_Lateral</i></p> <p><b>Define Orientation Using:</b><br/><b>Starting Point:</b> <i>Lab_Origin</i><br/><b>Existing Segment:</b> <i>Pelvis</i></p>                                                   | <ul style="list-style-type: none"> <li>Offset Using the Following <b>ML/AP/AXIAL</b> Offsets:<br/><b>X:</b> 0.1   <b>Y:</b> 0.0   <b>Z:</b> 0.0</li> <li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li> <li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li> </ul> |   |
| <p><b>2. Create Pelvis_Lateral_Project:</b></p> <ul style="list-style-type: none"> <li>Click <b>Landmarks</b> button</li> <li>Click <b>Add New Landmark</b> button</li> <li>Create Landmark: <i>Pelvis_Lateral_Project</i></li> </ul> | <p><b>Landmark Name:</b><br/><i>Pelvis_Lateral_Project</i></p> <p><b>Define Orientation Using:</b><br/><b>Starting Point:</b> <i>Lab_Origin</i><br/><b>Ending Point:</b> <i>Lab_X</i><br/><b>Project From:</b> <i>Pelvis_Lateral</i></p> | <ul style="list-style-type: none"> <li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li> <li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li> </ul>                                                                                                                   |  |

|                                                                                                                                                                                                                       |                                                                                                                                                                                                                     |                                                                                                                                                                                                                                                        |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <div><p><b>3. Create Lab_Lateral:</b></p><ul style="list-style-type: none"><li>Click <b>Landmarks</b> button</li><li>Click <b>Add New Landmark</b> button</li><li>Create Landmark: <i>Lab_Lateral</i></li></ul></div> | <div><p><b>Landmark Name:</b><br/><i>Lab_Lateral</i></p><p><b>Define Orientation Using:</b><br/><b>Starting Point:</b><br/><i>Lab_Origin</i><br/><b>Ending Point:</b><br/><i>Pelvis_Lateral_Projected</i></p></div> | <div><ul style="list-style-type: none"><li>Offset Using the Following <b>AXIAL</b> Offset: <i>0.1</i></li><li><b>DO NOT</b> Check: <i>Offset by Percent (1.0 = 100%)</i></li><li><b>DO NOT</b> Check: <i>Calibration Only Landmark</i></li></ul></div> | <div><div>Workspace   Signals and Events   Models   Reports   Cal</div><div>Segments   Landmarks   Muscles   Subject Data / Metrics</div><div>Digitizing   Pelvis Lateral   Pelvis Lateral Projected   Lab Lateral</div><div>Landmark Name: Lab_Lateral</div><div>Define Orientation Using:<br/>Starting Point Lab_Origin (Reference)<br/>Targets and/or Landmark:<br/>Ending Point Pelvis_Lateral_Projected<br/>*Lateral object<br/>*Project From<br/>* = Optional<br/>Existing Segment LAB</div><div>Landmark Offset from Start Point (Reference) or Segment Origin<br/>Offset to Existing Calibration Target or Landmark<br/>Offset Using the Following ML/AP/AXIAL Offsets<br/>ML<br/>AP<br/>AXIAL 0.1<br/>Offset by Percent (1.0 = 100%) (Meters when not checked)<br/>Calibration Only Landmark (Not generated for assigned motion file(s))</div><div>Undo Changes   Apply   Build Model   Close Tab</div></div> |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

Create Virtual Laboratory Segment

To create the Virtual Laboratory segment, switch to the **Segment** Tab in model builder mode and enter the following:

|                                                                                                                                                                                                                        |                                                                                                                                                                                                                                                                                                              |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <div><p><b>1. Create the Virtual Lab Segment:</b></p><p>1. In the Segment Name combo box type <b>Virtual Lab</b></p><p>2. Check the <b>Kinematic Only</b> Check Box</p><p>3. Select the <b>Create</b> button</p></div> | <div><div>Workspace   Signals and Events   Models   Reports   Cal</div><div>Segments   Landmarks   Muscles   Subject Data / Metrics</div><div>Segments   Segment Properties   IK Constraints</div><div>Segment Name<br/>Virtual Lab<br/>Segment Type<br/>Visual 3D<br/>Kinematic Only<br/>Create</div></div> |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

|                                                                  |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                       |                                                                                                     |                                                                                                |
|------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------|
| <div>4. In the <b>Virtual Lab</b> tab, enter these values:</div> | <div><b>Define Proximal Joint and Radius</b><br/><br/><b>Lateral:</b> <i>None</i>    <b>Joint:</b> <i>Lab_Z</i>    <b>Medial:</b> <i>None</i><br/><br/><b>Radius:</b> <i>0.01</i><br/><br/><b>Define Distal Joint and Radius</b><br/><br/><b>Lateral:</b> <i>None</i>    <b>Joint:</b> <i>Lab_Origin</i>    <b>Medial:</b> <i>None</i><br/><br/><b>Radius:</b> <i>0.01</i><br/><br/><b>Extra Target to define Orientation</b><br/><br/><b>Select:</b> <i>Lateral    Lab_Lateral</i><br/><br/><b>Select Tracking Targets:</b><br/><br/><b>Check:</b> <i>Use Calibration Targets for Tracking</i></div> | <div>5. Click on <b>Build Model</b>.<br/><br/>6. Click on <b>Close Tab</b> before proceeding.</div> | <div></div> |
|------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------|

We have just created a virtual lab that is aligned with the direction of progression of the subject.

Example: Virtual Lab Not Aligned with the Lab Axes

What if the subject is not walking along the lab axes? What if the subject is walking at a diagonal?

To calculate certain measures, such as the trunk angle relative to the lab, the virtual lab should be aligned with the subject's direction of progression. Although there are many methods to do this, some examples can be found [here](#).

From:  
<https://has-motion.com/wiki/> - Wiki Documentation Site

Permanent link:  
[https://has-motion.com/wiki/doku.php?id=visual3d:tutorials:modeling:virtual\\_laboratory&rev=1732899627](https://has-motion.com/wiki/doku.php?id=visual3d:tutorials:modeling:virtual_laboratory&rev=1732899627)

Last update: 2024/11/29 17:00

