

Language: English • français • italiano • português • español ****

==== Contents =====

- * [1 Introduction](#)
- * [1.1 Download](#)
- * [2 Prepare the Workspace](#)
- * [2.1 View files using the CMZ Library](#)
- * [3 Tag Files](#)
- * [4 Graph Data](#)
- * [4.1 1. Change Report Template to Plot Sub_01 tag](#)
- * [4.2 2. Overlay other subject's data](#)

Introduction

This tutorial is going to describe comparing different data sets. In this example, we're going to compare data between various subjects.

Download

The sample files can be downloaded [here](#).

- Data:
 - Sub 01 CMZ
 - Sub 02 CMZ
 - Sub 03 CMZ
- Scripts:
 - Tag.v3s (Example script that can be used to tag the files)
 - Tag_Loop.v3s (Example script that can be used to tag the files using a For_Each loop)
 - Overlay.v3s (Example script that can be used to overlay the report template with another subject)
 - Overlay_Loop.v3s (Example script that can be used to overlay the report template with another subject using a For_Each loop)
- Report Templates:
 - Report_Template_ALL_FILES.rgt
 - Report_Template_Sub_01.rgt

Prepare the Workspace


The first step is to load the data. It's recommended that you use the CMO Library to view multiple CMZ files at once. Many users simply choose to do a File → Open if they're viewing two small files. This option works for smaller data sets, but it's really recommended to use the CMO Library instead.

View files using the CMZ Library

1. **CMO Library → Set CMO Library Path**
 1. Select **..\Comparing Data Sets\Data folder**
 2. Click **OK**
2. **CMO Library → Browse CMO Library**
 1. Select **Yes**

[More detailed instructions to load the CMO Library can be found here.](#)

Tag Files

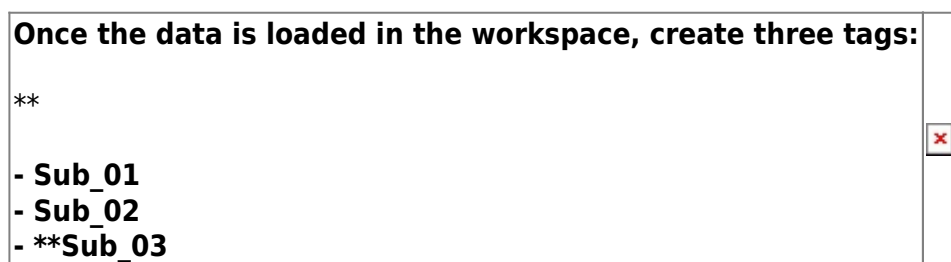

There are a few ways to create tags: (Choose one method of the following)

1. [Manually add tags](#) to the files or use the pipeline.
2. Use the [Assign_Tags_To_File](#) command to assign tags to subjects based on the filename. This can be done using the provided **Tag.v3s file**.
3. Using the [For_Each](#) command, you can loop over all the subjects in your workspace to create individual tags for them. This can be done using the provided **Tag_Loop.v3s file**.

Graph Data

Once the data has been loaded and tagged, create a report template that plots the data.

1. Change Report Template to Plot Sub_01 tag

A. **Open a report template** that plots the signals you are interested in viewing (ex. Open provided **Report_Template_ALL_FILES.rgt**) B. Change the tag from ALL_FILES (or whatever you're currently plotting) to one of the tags you want to plot (ex. Sub_01) and set the line color to the same for both the right and left sides. a. **Search/Replace:**

a1. Old File Tag: ALL_FILES	
a2. New File Tag: Sub_01	
a3. Click “Add»”	
a4. Old Line Color: #109900	
a5. New Line Color: Purple	
a6. Click “Add»”	
a7. Old Line Color: #990000	
a8. New Line Color: Purple	
a9. Click “Add»”	
a10. Select “Entire Report”	
a11. Click “Replace”	

D. **Save report template** (ex. Report_Template_Sub_01.rgt)

2. Overlay other subject's data


The step to overlay the subject data can be done manually (as described below), or using the pipeline (ex. provided scripts from download Overlay.v3s or Overlay_Loop.v3s).

A. Overlay Sub_02 on top of Sub_01 a. **Overlay Report:**

a1. Page No: ALL_PAGES	
a2. Identifier Type: TAG_NAME	
a3. Identifier Value: Sub_01	
a4. New Tag: Sub_02	
a5. New Color: Blue	

B. Overlay Sub_03 on top of Sub_01 a. **Overlay Report:** a1. Page No: **ALL_PAGES** a2. Identifier Type: **TAG_NAME** a3. Identifier Value: **Sub_01** a4. New Tag: **Sub_03** a5. New Color: **Cyan** C. Display the tag name in the Legend of the report

Retrieved from ""



From:
<https://has-motion.com/wiki/> - **Wiki Documentation Site**

Permanent link:
https://has-motion.com/wiki/doku.php?id=visual3d:tutorials:reports:comparing_data_sets&rev=1718644148

Last update: **2024/06/17 17:09**

